

KINOS X

Focus on your patients



LOJER®

For Easy
Care





Kinos X5



Kinos X4

Main features:

1. Can be used in prone, supine and sitting positions
2. Sturdy construction that can handle even the heaviest patients
3. Foot bar adjustment leaves hands free

Lojer Kinos X treatment table is designed for professionals working in the demanding and diverse rehab world. The Kinos X has a wide height adjustment range, making it easy to adjust for different needs. The table is now equipped with electrically raised castors as a new feature, making it incredibly easy to move around. Lojer has done a lot of testing and designing so that you can work stress-free. The Kinos X Treatment table has passed all the medical device regulation tests performed by an external testing laboratory.

The range:

- Kinos X5 has a 65 cm wide two-piece table-top with back section adjustment. The back section and height of the table are hygienically adjusted with foot bars on both sides of the table.
- Kinos X4 has a 55 cm wide one-piece table-top with side supports. The height of the table is hygienically adjusted with a foot bar on both sides of the table.



The back section and height of the table are hygienically adjusted with a foot bar, leaving your hands free to help and support the patient.

Lojer Kinos X treatment table is sturdy and strong with High Safe Working Load.




The upholstery material, Future, is antimicrobial by ISO-certification. The material is easy to clean and it's six time more durable than artificial leather.



Technical data and features

Model	X4	X5
Total width	55 cm	65 cm
Total length	203 cm	203 cm
Total weight	98 kg	109 kg
Height adjustment	46 - 95 cm	46 - 95 cm
Back section length	-	52 cm
Back section adjustment	-	0 - 70°, electrical
Head section adjustment	-90°...+35° (U-model)	-90°...+35° (U-model)
Castors	Ø65 mm, electric lifting	Ø65 mm, electric lifting
Other standard features	Safety feature: auto-lifting in emergency situations Safety feature: table will be automatically locked after 120 minutes if not in use Safety feature: motors do not pull, but lower using gravity Arm rests Side supports One long height adjustment bar on both sides of the table Electric lifting on castors	Safety feature: auto-lifting in emergency situations Safety feature: table will be automatically locked after 120 minutes if not in use Safety feature: motors do not pull, but lower using gravity Arm rests Two height adjustment bars on both sides of the table for height and back section Electric lifting on castors Moulded plug for head section
Upholstery	Material Future, color Carbon	Material Future, color Carbon
Frame color	Carbon	Carbon
Safe working load (SWL)	210 kg	210 kg
Electric system	24V, IPX6	24V, IPX6
CE marked	Yes, MDR	Yes, MDR
Warranty	2 years	2 years
Country of Manufacture	Finland	Finland

Kinos X4

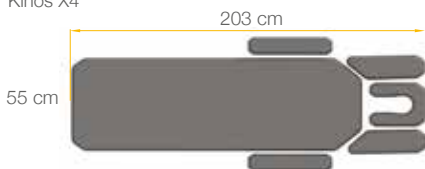

Kinos X5

